

## LITHOSTAR Multiline

**SP**

### Installation and Setting Instructions

DICOM Connect

© Siemens AG 2003

The reproduction, transmission or use of this document or its contents is not permitted without express written authority. Offenders will be liable for damages. All rights, including rights created by patent grant or registration of a utility model or design, are reserved.

English

Doc.Gen.Date: 02.03

Print No.: RXL2-120.814.02.01.02

Replaces: n.a.

**Document revision level**

The document corresponds to the version/revision level effective at the time of system delivery. Revisions to hardcopy documentation are not automatically distributed.

Please contact your local Siemens office to order current revision levels.

**Disclaimer**

The installation and service of equipment described herein is to be performed by qualified personnel who are employed by Siemens or one of its affiliates or who are otherwise authorized by Siemens or one of its affiliates to provide such services.

Assemblers and other persons who are not employed by or otherwise directly affiliated with or authorized by Siemens or one of its affiliates are directed to contact one of the local offices of Siemens or one of its affiliates before attempting installation or service procedures.

**Application to Systems**

DICOM can be built into the LITHOSTAR Multiline only if no MULTISPOT is installed.

**Required Tools**

- Tool case (see "Service Tools Catalogue" in the "Spare Parts Catalogue").
- No special tools are needed.

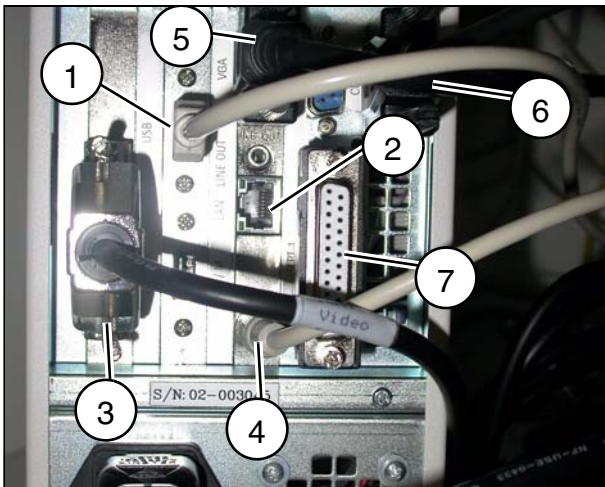


Fig. 1

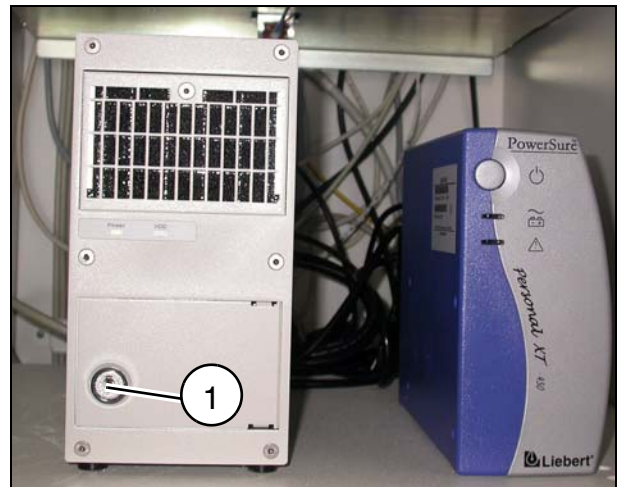


Fig. 2

## Procedure

- Open the control console door.
- Place the monitor, keyboard and CD drive in a suitable place on the control console.
- Place DICOM and UPS (uninterrupted power supply) below the control PC in the console (Fig. 2).
- Run the keyboard, USB and monitor extension cables to DICOM.
- Connect the CD writer to the USB (1/Fig. 1).
- Connect the network cable (2/Fig. 1).
- Connect the video cable (3/Fig. 1) and screw it tight. Connect the other end of the cable to the BAS cable designated "Multispot" via the supplied BNC connector.
- Connect the keyboard cable (4/Fig. 1).
- Connect the plug for the video signal monitor (5/Fig. 1) and screw it tight.
- Connect the connecting cable DICOM - UPS (6/Fig. 1) and screw it tight.
- Plug the dongle (7/Fig. 1) into the LPT1/Printer socket and screw it tight.
- Run the mains cable from the UPS (uninterrupted power supply) to the multipoint connector in the console and connect it.
- Connect the power supply for CD drive to the power strip in the console.
- Connect the DICOM and monitor mains cables to the UPS sockets.
- Open the front panel (1/Fig. 2) with the key and check whether the power switch underneath it is switched on.



Fig. 3



Fig. 4

- Plug the cables coming from the keyboard as shown in figure 3 and connect them to the keyboard extension cable.
- Run the control console trigger cable designated "Hardcopy" to the keyboard and connect it to the trigger cable of the keyboard as shown in figure 4.

### DICOM start-up

Start-up is carried out as described in the supplied instructions.